SEFAR® CleanDisc

Inspection of Technical Cleanliness






SEFAR® CleanDisc System – Mesh filter

Modern laboratories are using mesh filters to perform the Inspection of Technical Cleanliness according VDA 19-1 & ISO 16232.




Extraction and Filtration

- Good chemical compatibility when using common extraction liquids.
- Increase of throughput rate and defined particle retention rate, in comparison to membranes.

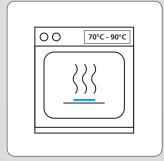


Drying

- Improved flatness and reduced rolling behavior by selecting the appropriate polymer type (PET recommended) adaption of drying time and temperature.
- PET: max. 30 min@ 90°C/194°F
- PA: max. 30 min@ 70°C/158°F

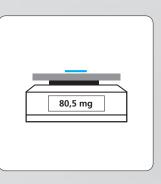

Gravimetrical Analysis

- Polyester offers less humidity absorption behavior compared to Nylon (PA) or cellulose membranes.
- Performing gravimetrical analysis, Polyester is of advantage, as it is less sensitive to temperature or humidity fluctuations.

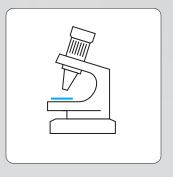

Light optical analysis

- Mesh filters have a lower filter occupancy compared to membranes, which improves the analyzability.
- Yellow mesh filters allow the detection of white particles and increase the contrast of dark particles.





SEFAR



used for Inspection of Technical Cleanliness

Clean, safeguarded and easy to handle mesh filter, dispensed to directly perform a technical cleanliness analysis.



SEFAR® CleanDisc Dispenser

- Easy and single withdrawal of the filter disc.
- Enclosed design avoids cross contamination from the environment.
- Storage in a safeguarded environment, according to the norm.



SEFAR® CleanDisc Cassette

- Easy exchangeable cassettes allow operating one dispenser with various cassettes.
- Clear identification of the operated mesh filter type displayed on the label e.g. 5 μm PET mesh filter.



SEFAR® CleanDisc mesh filter

- The safety foil secures the level of cleanliness until usage.
- Mesh filter openings correspond to the defined particle size classes in the norm.
- Laser cut edges avoid frying edges and loose yarn threads.
- Produced in cleanroom conditions.




SEFAR® CleanDisc order options

- Filter mesh opening [μm]: 1, 5, 8, 10, 15, 20, 40, 50, 80, 100, 150
- Diameter: 47 mm or 25 mm
- Material: PET or PA (Nylon)
- Color: white (yellow on request)
- Package size: 200 or 600 pieces



Developed with experts,

tested by users.

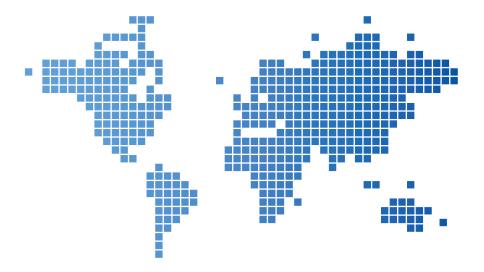

Sefar Worldwide

Sefar is the leading manufacturer of monofilament precision fabrics for filtration and screen printing applications. Sefar products are used in a wide variety of industries, reaching from electronics, medical, automotive, aerospace, food, water, pharmaceutical and graphics to chemicals, fertilizers and minerals. With its profound understanding of the applications, Sefar helps its customers to achieve optimum results in their industrial processes. With subsidiaries and representatives around the world, Sefar offers its customers the comprehensive, global support they need.

Sefar AG

Hinterbissaustrasse 12 9410 Heiden Switzerland

Phone +41 71 898 57 00

CleanDisc@sefar.com www.cleandisc.com

www.sefar.com

Africa Contacts

South Africa – Johannesburg Phone +27 11 708 2485 info@sefarfilterpure.com

Morocco – Casablanca Phone +212 522 66 14 66 info.ma@sefar.com

Americas Contacts

USA – Buffalo Phone +1 800 995 0531 sales@sefar.us

Canada – Saguenay Phone +1 418 690 0888 info.bdh@sefar.com

Brazil – São Bernardo do Campo Phone +55 11 4390 6300 contato.brasil@sefar.com

Mexico – Cuautitlan Izcalli Phone +52 55 2291 0652 info.mexico@sefar.com

Oceania Contacts

Australia – Huntingwood Phone +61 2 8822 1700 oceaniasales@sefar.com

New Zealand – MT Wellington Phone +64 9 527 4409 salesnz@sefar.com Asia Contacts

China – Suzhou Phone +86 512 6283 6383 info.cn@sefar.com

India – District Thane Phone + 91 9028 656 256 info.in@sefar.com

Malaysia – Johor Phone +60 7 861 0023 info.my@sefar.com

Singapore

Phone +65 6299 9092 info.sg@sefar.com

South Korea – Seoul Phone +82 2 713 9133 info.kr@sefar.com

Thailand – Bangkok Phone +66 2 210 0018 info.th@sefar.com

Vietnam – Ho Chi Minh City Phone +84 862 925 562 info.vn@sefar.com **Europe Contacts**

England – Bury Phone +44 161 765 3530 sales@sefar.co.uk

France – Billy-Berclau Phone +33 3 21 08 70 70 vente@sefar.com

Germany – Edling Phone +49 8071 90400 0 info@sefar.de

Italy – Collegno Phone +39 011 42 001 info@sefar.it

Netherlands – Lochem Phone +31 573 744 111 info@sefar.nl

Poland – Poznan Phone +48 511 467 467 biuro.pl@sefar.com

Spain – Cardedeu Phone +34 93 844 47 10 info.maissa@sefar.com

Turkey – Ankara Phone +90 312 287 34 14 info.tr@sefar.com

